
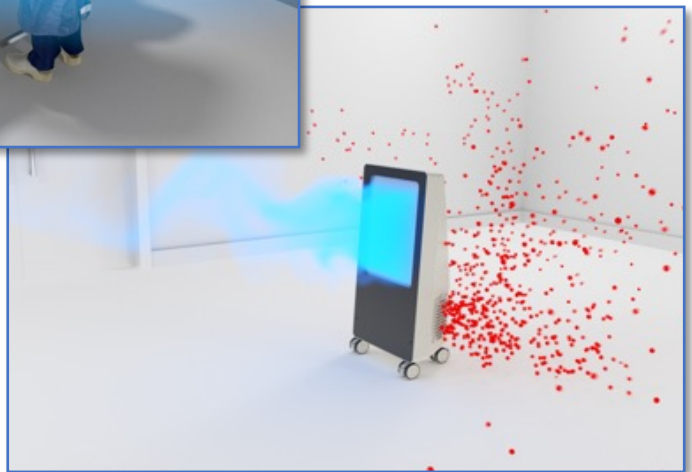


Operio CLINIC

Ventilation system for veterinary clinic.

Mobile ultra-clean air production unit for the ventilation system of veterinary clinics.

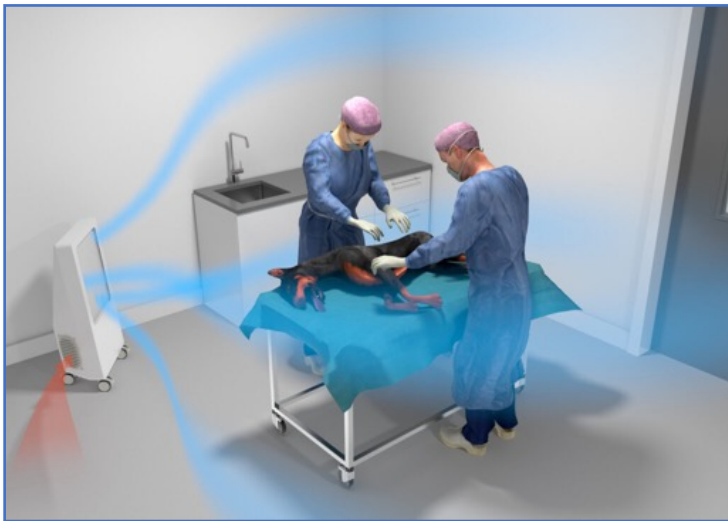

TOUL Operio CLINIC

- Eliminates 99.995% of viruses, bacteria, fine dust, etc.
- Controlled mechanical ventilation for veterinary clinics.
- Create overpressure (positive pressure) in the room.
- Increase the number of air changes/hour.

Normeditec

What is Operio CLINIC?

Operio CLINIC is a mobile air purification unit intended for use in any indoor environment. It produces a continuous and ultra-clean air flow.

Operio CLINIC is a mobile device for air disinfection based on HEPA 14 filters which eliminates 99.995% of all airborne, bacterial and viral particles, dust, etc.

The filter also has an internal UVC function, which effectively kills all filtered bacteria and viruses, including the corona virus, whenever we want.

The Operio CLINIC device can be used to clean ambient air in an enclosed space or directed towards a specific area. No system work is necessary for installation: the device only requires a common power socket.

TECHNICAL SPECIFICATIONS	
Dimensions	11 x 47 x 97 cm
Weight	15 Kg
Air cleaning capacity / hour: selection of different fan speeds	Speed 1: 200 m ³ /h Speed 2: 300 m ³ /h Speed 3: 480 m³/h
HEPA 14 filter	Filters 99.995% particles diameter >0.1 micron
UV light	Manual activation with automatic switched off after 1 min.



Normeditec

Clinical applications

Ventilation for veterinary clinics and operating environments.

Operio CLINIC can be used to greatly improve aseptic conditions and perform minor veterinary surgery in small outpatient settings.

Opero CLINIC can be used to improve the level of asepsis during veterinary procedures in operating room and elsewhere. The device has different speeds and can be used to achieve the desired degree of asepsis for different types of surgery. Installation takes only a few seconds and can be transported between operating rooms and other environments.

**Operio CLINIC:
High aseptic level with limited budget.**



Overpressure and number of air exchange/hour.


Controlled mechanical ventilation for veterinary clinics.

The device sanitizes 480 m³ of air every hour.
For a standard room of 25m² (2.7 meters high) this means
the device filters the air in the room every 7 minutes.

CREATE POSITIVE PRESSURE and INCREASE the number of SPARE PARTS/HOUR.

Operio CLINIC is a mobile device for air sanitization based on HEPA 14 filters which eliminates 99.995% of bacterial, viral particles, dust, dust and any airborne particle. Operio Clinic can be connected to external air to provide fresh air into the room and to create excess pressure compared to the surrounding environments.

The simplest and most economical way to create overpressure or increase the number of air changes.

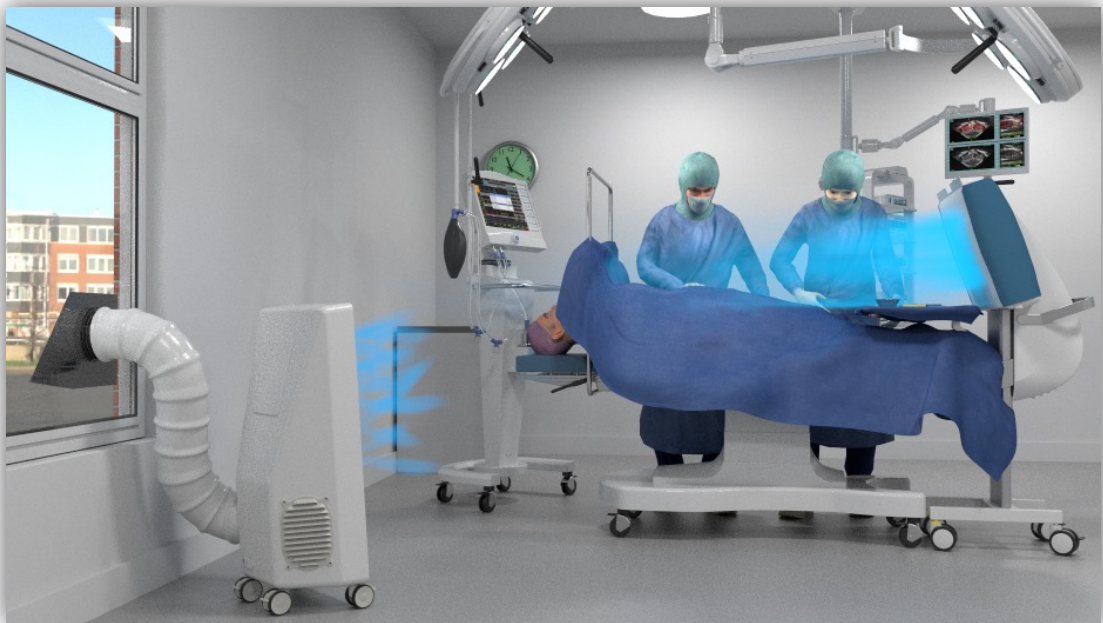

NORMEDITEC S.R.L.

Via De Gasperi 19 43018, Sissa Trecasali

Tel.0521 878 949 Tel.348 730 2445 – 331 185 9110 Fax.0521 363 731

www.normeditec.com info@normeditec.com

